Primado Air







To us, KANUMA is the center of the world.

It was here in Kanuma that we first started and it is here where we will continue to create the best in class, high-speed rotation technology in the world. "Made in Japan" symbolizes world class quality and craftsmanship and now "Made in Kanuma" will define the finest quality in high-speed drills across the world.





Quality Results from Processes

This has been the NSK policy since its foundation. Up to 30 processes may be necessary to create one complete component from a bar of raw material. The total number of machining processes required to create all parts of a single product is 930. If the individual processes do not attain the degree of accuracy needed, the final product will not deliver the requisite high precision. Ensuring micron-scale accuracy is crucial. In pursuit of such accuracy, NSK's key strength lies in the ability to internally manage integrated production from the component stage. To ensure high quality, defective parts are never passed. Manufacturing is a constant battle to push the limits. And the state-of-the-art manufacturing pursued by NSK is an never-ending challenge.



Global Network

NSK has sales subsidiaries in 15 countries with each subsidiary actively conducting sales and marketing activities in its local territory, providing detailed after-sales services and attracting new NSK fans every day. Furthermore, NSK is also collaborating with import agents in other countries and building strong ties to promote sales in local areas. NSK has built a solid global sales network and exports to 135 countries. The voice of our sales partners is important, and NSK strives to ensure quick and reliable support. These strong partnerships can then transfer to our customers, who can be fully confident in our brand and products.



Stablity Of Silence

Stable Rotation


- •The special pneumatic motor provides both powerful cutting performance and low noise
- ·Less run-out delivers precise, smooth bone cutting performance

Uncompromising Design

- · Variable exposure allows for versatility of surgical procedures *P300 Attachment
- ·High quality cutting accessories made with innovative technology


Better Usability


- •The universal bur length system features one bur shank length that can be used with any length attachment *standard Attachment 300
- •The risk of excessive heat development and consequent malfunction as a result of insufficient cleaning is reduced

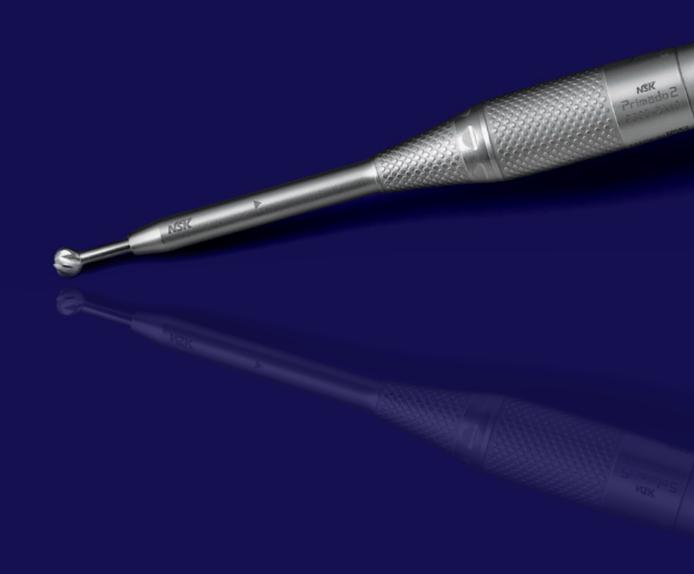



Product Lineup

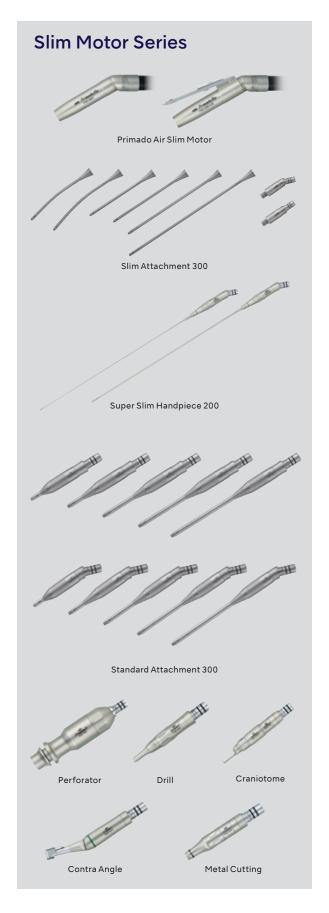

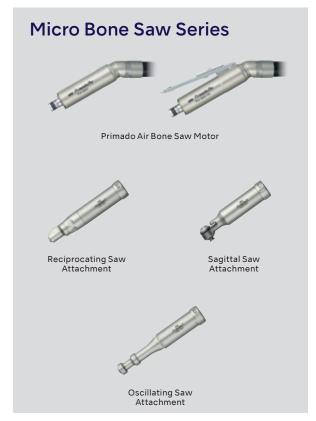





Primado Air Foot Control

The incorporated relief valve protects against overpressure and prevents excessive speed

An open pedal design enables precise and easy control

A control mode switch lever allows either foot or hand operation modes



Order Code
PDA-FC
Primado Air Foot Control

• Variable speed control

PDA-FCH-SCH
Foot Control Hose/Schrader



PDA-FCH-DN
Foot Control Hose/N-DISS



Slim Motor Series

Compact and lightweight design prevents surgeon fatigue during procedures

A variety of attachments are available

Quick Joint allows for an easy attachment-motor connection



Primado Air Slim Motor

- Speed: Max. 80,000 min⁻¹
 Bundled Item: PDA-CP-SMH, PDA-CP-MP (see p.19)
 PDA-EZ-SMH, P200-EZ-SM (see p.20)



PDA-SMH Primado Air Slim Motor

PDA-SMH-HS

Primado Air Slim Motor/HS · Hand Control type


P300 Attachment Key Features



The stable rotation performance is consistent with current P200 attachments, however, the P300 attachments are designed to improve visibility, smooth operation, and safety.

Variable Exposure

• Improved accessibility and visibility due to the 6-step adjustment mechanism

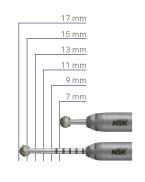



25.5 mm

23.0 mm

20.5 mm

18.0 mm

15.5 mm

13.<u>0 mm</u>


Slim Attachment 300

 $\bullet \ \, \text{The attachment's twist collet allows for variable exposure without having to handle the bur.} \ ^* \text{Standard Attachment 300} \\$



Universal Bur Length System *Standard Attachment 300

- One bur size fits all lengths of attachment
- Easier management in the operating room

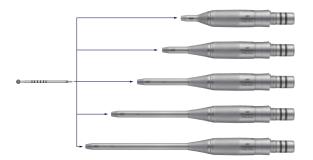


Slim Attachment 300 80,000min⁻¹

Slim Attachment Hub






Slim Tube

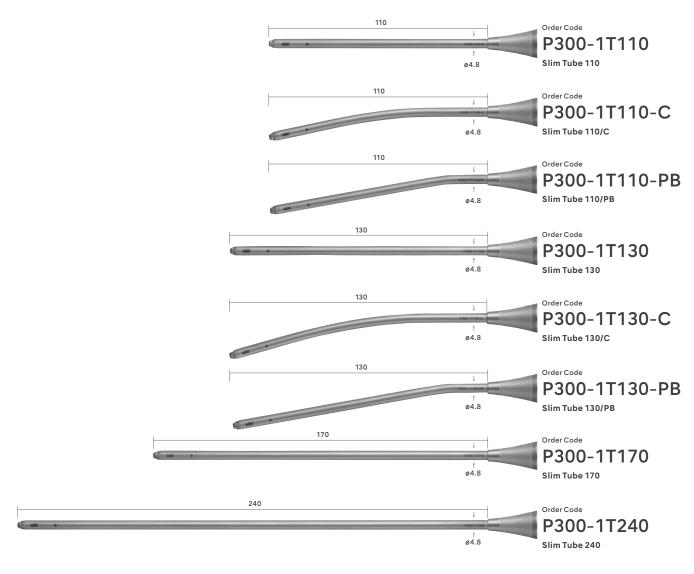



Standard Attachment 300 80,000min⁻¹

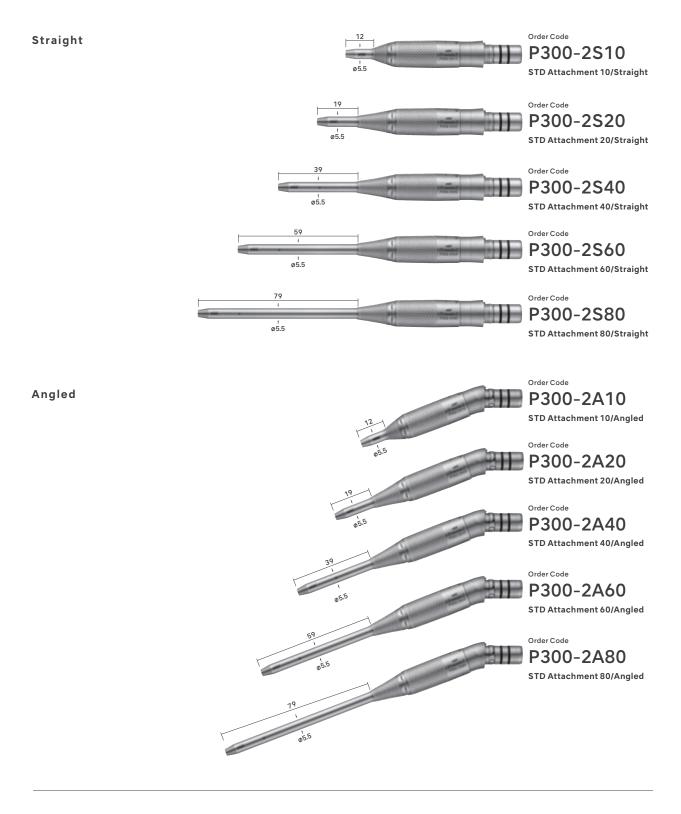




P300-2-B
STD Attachment Beak




Slim Handpiece 200 80,000min⁻¹



Super Slim Handpiece 200 20,000min⁻¹

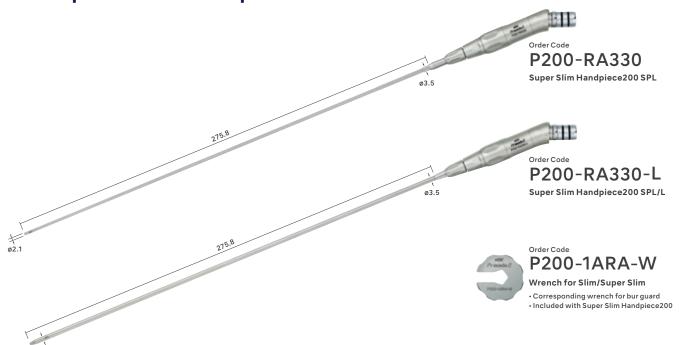

Standard Handpiece 200 80,000min⁻¹

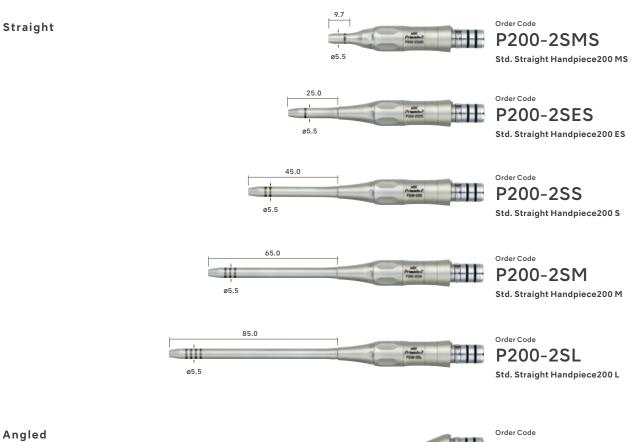





Perforator Handpiece 200 1,250min⁻¹



Order Code P200-PER Perforator Handpiece200
• Compatible Chuck: Hudson Chuck

Drill Handpiece 200 80,000min⁻¹



Order Code P200-2SD Straight Drill Handpiece 200

Contra Angle Handpiece 200



Order Code

P200-CA-4

Contra Angle Handpiece 200 4:1

- 4:1 Reduction

- (Speed: Max. 20,000 min⁻¹)
 Latch Type Chuck
 Corresponds to Ø2.35mm CA burs

Craniotome Handpiece 200 80,000min⁻¹



P200-CRA
Craniotome Handopiece200



P200-CRD-P
Pediatric Duraguard 200



P200-RCRD-P
Rotatable Duraguard200 P



P200-CRD-M
Medium Duraguard200



P200-RCRD-M
Rotatable Duraguard200 M



P200-CRD-L
Long Duraguard200



Order Code
P200-RCRD-L
Rotatable Duraguard200 L

Metal Cutting Handpiece 200 80,000min⁻¹



P200-3MC
Metal Cutting Handpiece200

Micro Bone Saw Series

A modular system in which one motor accepts 3 different types of saws

Designed for powerful cutting performance and less vibration



Primado Air Bone Saw Motor

- Speed: Max. 20,000 min⁻¹
 Bundled Item: PDA-CP-BMH, PDA-CP-MP (see p.19)
 PDA-EZ-BMH, P200-EZ-BM (see p.20)



Order Code PDA-BMH

Primado Air Bone Saw Motor

Order Code PDA-BMH-HS Primado Air Bone Saw Motor/HS · Hand Control type

Reciprocating Saw 20,000min⁻¹



Sagittal Saw 20,000min⁻¹



Oscillating Saw 20,000min⁻¹



Accessories

Cleaning Tool

The risk of excessive heat development and consequent malfunction as a result of insufficient cleaning is reduced





Manual cleaning

Order Code
PD-CTH
Cleaning Holder

• Bundled Item: PD-CTTU-606501
PD-CTTU-600504



Order Code
PD-CTAD-SM
Cleaning Adaptor SM
For Slim Motor Series each attachment



PD-CTAD-1T
Cleaning Adaptor 1T
For Slim Tube



Order Code
PD-CTTU-606501
Cleaning Tube 6.0/L650



Washer disinfector

Order Code
PD-CTAD-G
Cleaning Adaptor G
• For drying




Order Code
PD-CTAD-LM
Cleaning Adaptor Luer M



PD-CTAD-LF
Cleaning Adaptor Luer F



For automated cleaning For manual cleaning For drying

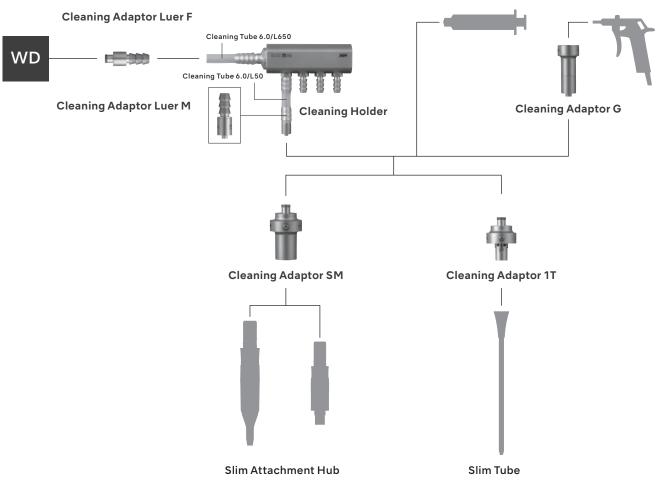

Slim Motor Series each attachment

^{*}Please be sure to read the reprocessing manual.

^{*}These Cleaning Tools can not be used with the Bone Saw Attachment.

Sterilization Case



Order Code SC-ML Sterilization Metal Lid



Order Code SC-T-F Sterilization Tray/Free • Dimensions : W480 x D250 x H60 (mm)



Order Code

SC-C-H70

Sterilization Case/H70

- Dimensions : W512 x D262 x H70 (mm) Single level case



Order Code

SC-SM

Silicone Mat

• The mat for the Sterilization Tray



SC-C-H125

Sterilization Case/H125

- Dimensions : W512 x D262 x H125 (mm)
 Double level case

Cap



Order Code

PDA-CP-SMH

Cap/Slim Motor



Order Code PDA-CP-MP Cap/Motor Plug



Order Code

PDA-CP-BMH

Cap/Bone Saw Motor

Lubrication





Order code


Z182100

PANA Spray Plus
• For attachments, 480 mL

Order code

Y900630

PANA Spray Plus (6pcs)
• Pack of 6

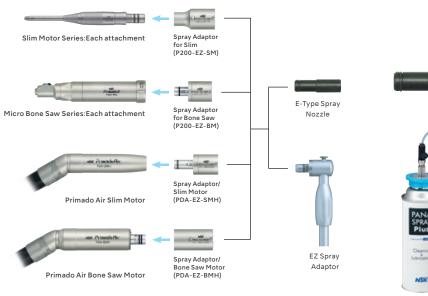

ASK

ASK

Order code

Z019090

E-Type Spray Nozzle

Order code

PD-EZ-M

EZ Spray Adaptor
• PANA Spray Plus is not included

Spray Mist Absorber





Order Code **Y900084**

Spray Mist Absorber

• Pack of 10

Index

Order Code	Page
P200-1A135	P.11
P200-1A155	P.11
P200-1A155-C	P.11
P200-1A175	P.11
P200-1A175-C	P.11
P200-1ARA-W	P.11
P200-2AES	P.12
P200-2AL	P.12
P200-2AM	P.12
P200-2AMS	P.12
P200-2AS	P.12
P200-2SD	P.13
P200-2SES	P.12
P200-2SL	P.12
P200-2SM	P.12
P200-2SMS	P.12
P200-2SS	P.12
P200-3MC	P.14
P200-CA-4	P.13
P200-CRA	P.14
P200-CRD-L	P.14
P200-CRD-M	P.14
P200-CRD-P	P.14
P200-EZ-BM	P.20
P200-EZ-SM	P.20
P200-OSC	P.16
P200-PER	P.13
P200-RA330	P.11
P200-RA330-L	P.11
P200-RCRD-L	P.14

Order Code	Page
P200-RCRD-M	P.14
P200-RCRD-P	P.14
P200-REC	P.16
P200-SAG	P.16
P300-1AHA	P.09
P300-1AHS	P.09
P300-1T110	P.09
P300-1T110-C	P.09
P300-1T110-PB	P.09
P300-1T130	P.09
P300-1T130-C	P.09
P300-1T130-PB	P.09
P300-1T170	P.09
P300-1T240	P.09
P300-1T-B	P.09
P300-1T-H	P.09
P300-2A10	P.10
P300-2A20	P.10
P300-2A40	P.10
P300-2A60	P.10
P300-2A80	P.10
P300-2-B	P.10
P300-2S10	P.10
P300-2S20	P.10
P300-2S40	P.10
P300-2S60	P.10
P300-2S80	P.10
PDA-BMH	P.15
PDA-BMH-HS	P.15
PDA-CP-BMH	P.19

Page
P.19
P.19
P.20
P.20
P.06
P.06
P.06
P.07
P.07
P.17
P.20
P.16
P.19
P.20
P.20
P.20
P.20



Create it